


GenXFlow™

Customer User Manual

Genomics Pipeline Management Platform

Version	1.0.0
Effective Date	August 2026
Platform	Windows Desktop
Audience	GenXFlow Customers

This manual explains how to install, activate, configure, run, monitor, and manage GenXFlow. It is written for customers and day-to-day users of the application.

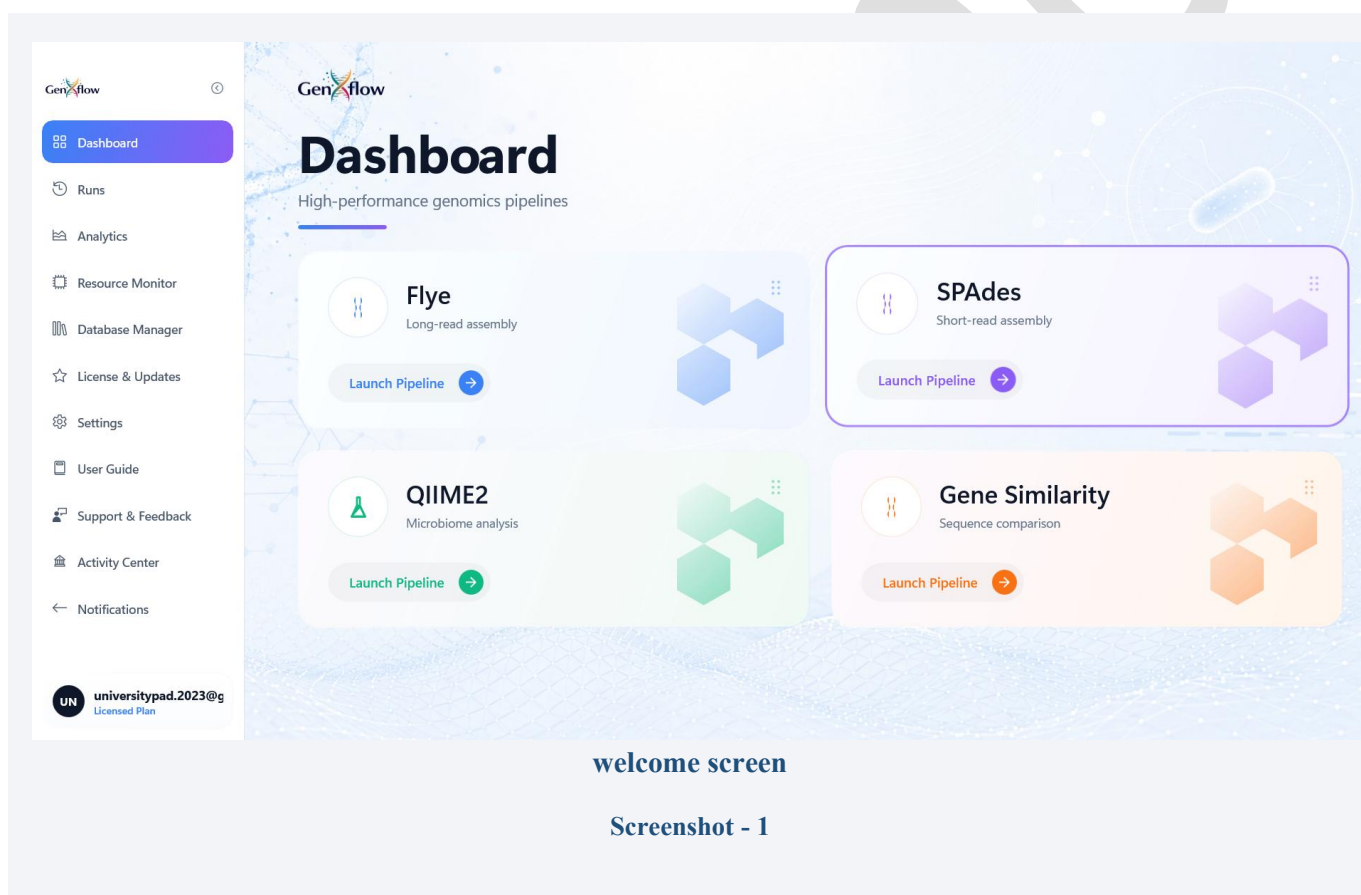

Contents

1. Before You Start
 2. Installation & First Launch
 3. Account & Sign-In
 4. License & Subscription
 5. Dashboard & Navigation
 6. Running Genomics Workflows
 7. Monitoring & Run Management
 8. Reference Databases
 9. Notifications & System Health
 10. Application Updates
 11. Feedback & Troubleshooting
 12. Support Notes
- Appendix A. Screenshot Checklist

1. Before You Start

GenXFlow is a Windows desktop application for managing common genomics analysis workflows, including genome assembly, amplicon analysis, sequence similarity analysis, and BLAST searches.

1.1 Recommended system requirements

Component	Minimum	Recommended
Operating system	Windows 10 64-bit	Windows 11 64-bit
Memory	8 GB RAM	32 GB RAM or more
CPU	4 cores	16+ cores
Disk space	50 GB free	500 GB+ SSD
Internet	Required for authentication, licensing, downloads, and updates	Stable broadband connection

1.2 Required software

Software	Purpose
Docker Desktop	Required for SPAdes, Flye, and QIIME2 workflows
GenXFlow	Main customer application; supplied by ArqGene


Before running a Docker-based workflow, make sure Docker Desktop is installed and available. GenXFlow can attempt to start Docker Desktop when needed.

2. Installation & First Launch

Use the installer supplied by ArqGene or your organization. Do not use an installer obtained from an unofficial source.


Make sure docker desktop is installed and opened before starting the installer


1. Run the GenXFlow installer.
2. Accept the license agreement.
3. Choose the installation directory, or keep the default location.
4. Select Install and wait for the installation to complete.
5. Select Finish. GenXFlow will open automatically if enabled by the installer.

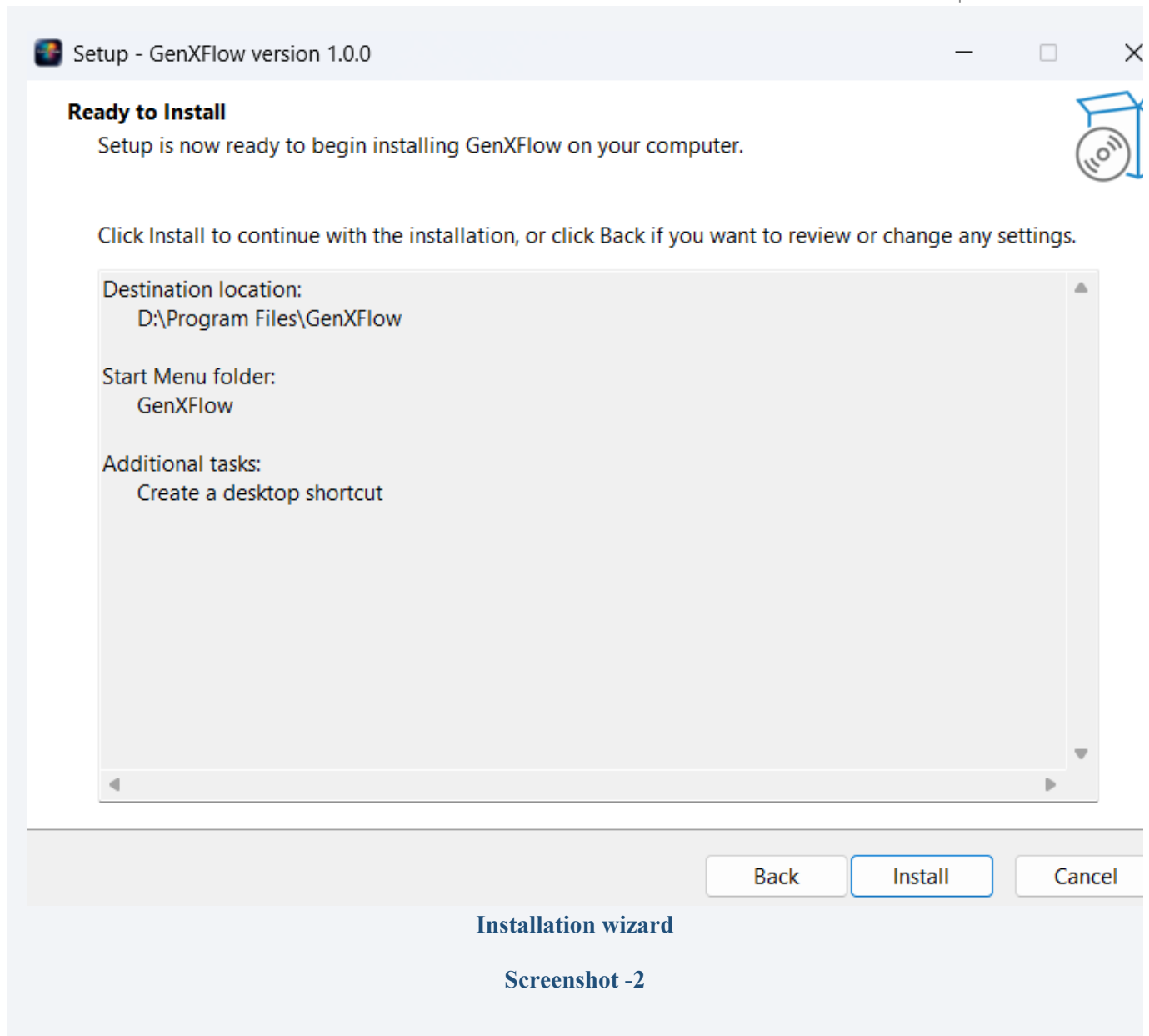

2.1 First launch

During the first launch, GenXFlow initializes the application, checks account/license status, checks Docker availability, and loads notifications or pending activities.

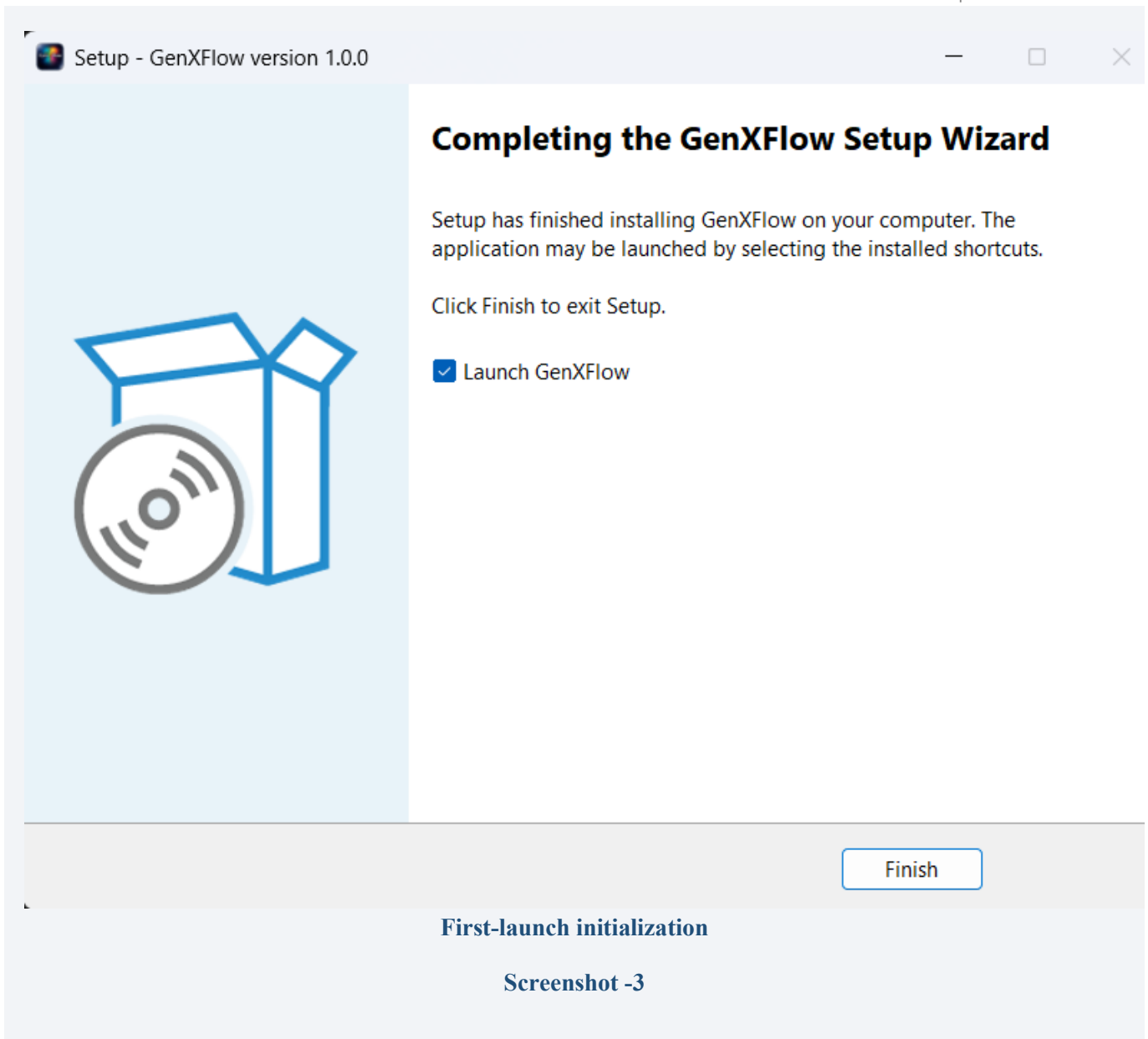

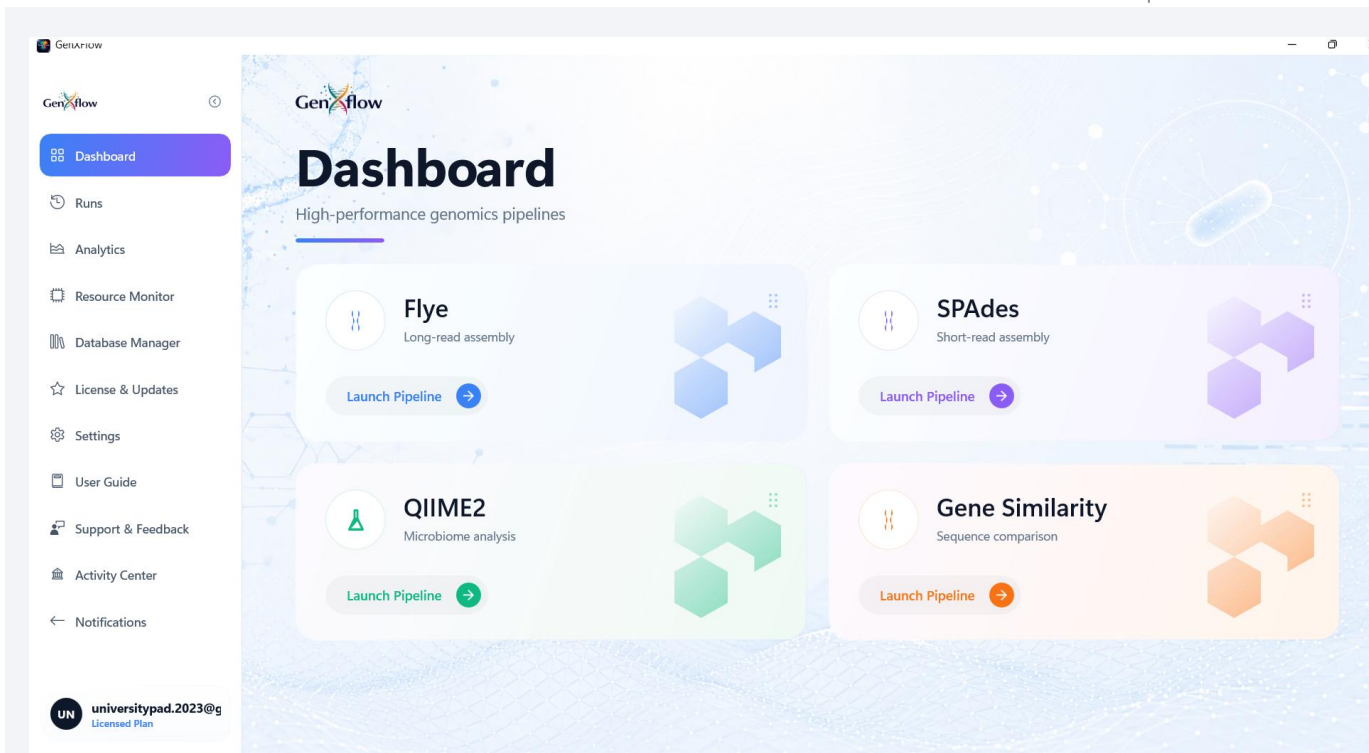


Screenshot -4

3. Account & Sign-In


3.1 Sign in

1. Open GenXFlow from the desktop shortcut or Start Menu.
2. Enter your registered email address and password.
3. Select Sign In.
4. After successful authentication and license validation, the main dashboard opens.

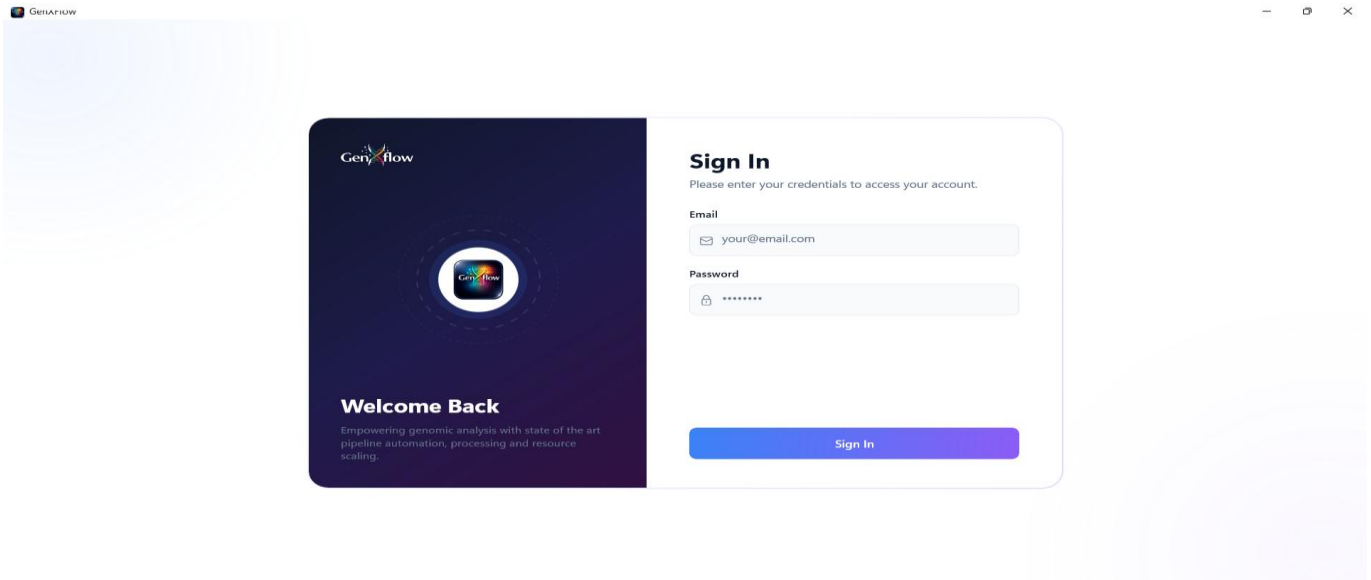

Login page

Screenshot - 5

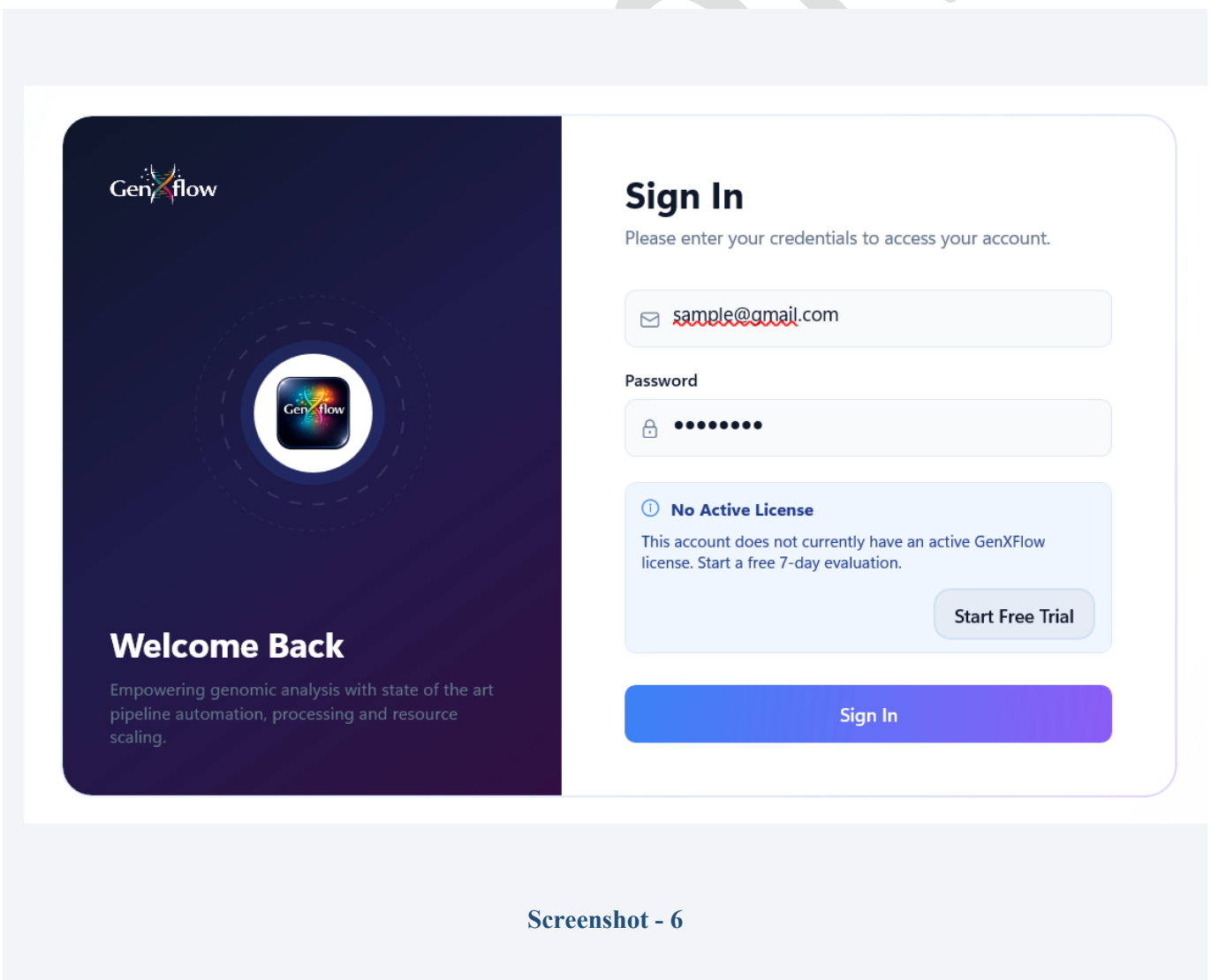

Screenshot - 6

3.2 Create an account

1. Select Create Account or Register from the sign-in screen.
2. Enter your username, email, password, and any optional recovery information.
3. Select Register.
4. Return to the sign-in screen and sign in with the new account.

3.3 Free trial

Where enabled by your customer package, the free trial provides 2 pipeline executions for 7 days. Once the trial allocation is exhausted, a paid license is required for further pipeline runs.

3.4 Password changes

For security reasons, GenXFlow does not provide a self-service Forgot Password or password-reset workflow in the customer application. If you need to change or reset your password, contact ArqGene support or the support channel supplied with your customer license. Do not send your existing password in a support request.

4. License & Subscription

4.1 License types

License	Description	Pipeline runs
Trial	Limited evaluation access	2 total for 7 days
Professional	Paid individual license	Unlimited
Enterprise	Paid organizational license	Unlimited

4.2 View license status

1. Open Account or Settings.
2. Open the License or Subscription section.
3. Review status, edition, registered email, expiry date, device information, and remaining trial runs where applicable.

License & Updates
Manage workstation activation, view module installations checklist, and pull latest offline system updates.

License Details Active

Edition:	GenXFlow Professional Edition
Licensed To:	universitypad.2023@gmail.com
User ID:	43f351d9-8f24-4e48-b5ee-f813441d0afd
License Type:	Licensed
License Expiration:	Aug 01, 2126 (36500 days remaining)
Device ID:	e822fcb1dacdded885defa83e691ac8d71282a97163abd74cbfdbc52b9f49a34

[Deactivate This Workstation](#)

Workstation Pipelines & Modules
Status of optional pipeline components locally configured in the app.

☒ Assembly Module Installed & Active	☒ QIIME2 Module Installed & Active	☒ FLYE Module Installed & Active
--	--	--

Software Version & Updates


Current Version: 1.0.0	Latest Version: 1.0.0	Up to Date
------------------------	-----------------------	------------

[Check for Updates](#)

SCREENSHOT 7
License information page

GenXFlow validates license eligibility when you start a pipeline. An internet connection may be required for account and license validation.

5. Dashboard & Navigation

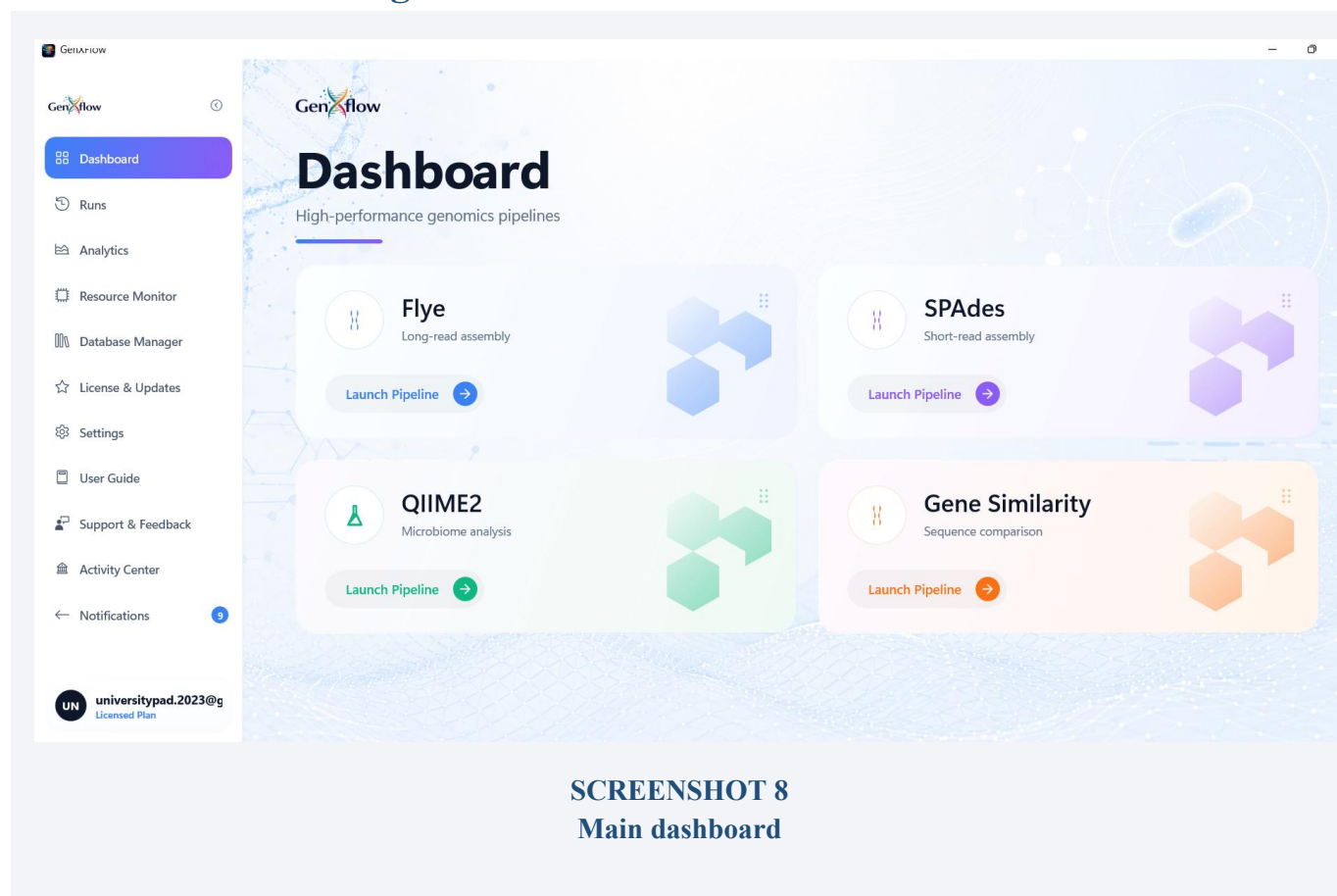

The dashboard provides quick access to the major workflow areas:

- Pipeline cards for SPAdes, Flye, QIIME2, Similarity Analysis, and BLAST.
- Run statistics and recent activity.
- System health information such as CPU, RAM, disk, and Docker status.
- Notifications and account/settings access.

5.1 Main navigation

Menu	Purpose
Dashboard	Overview of activity and system status
Pipelines	Open an analysis workflow
Run History	Review active and completed runs
Databases	Download and manage reference databases
Notifications	View system and pipeline notifications
Settings / License	Manage account, license, feedback, and application settings

6. Running Genomics Workflows

For SPAdes, Flye, and QIIME2, Docker Desktop must be available before execution. Select the workflow, provide the required input files, configure the available parameters, choose an output directory, and start the run.


6.1 SPAdes Genome Assembly

Use SPAdes for short-read and hybrid/long-read genome assembly workflows supported by your GenXFlow configuration.

Inputs and settings

Item	Description
Short R1 / R2	Paired short-read FASTQ inputs when using short-read data
Long Reads	Long-read FASTQ input for supported configurations
Output Directory	Location for generated results
SPAdes Type	Assembly mode such as standard, rna, plasmid, meta, or other supported modes
Genome Size	Estimated genome size, when required
Threads	Number of CPU threads to use
Pipeline Steps	Tools such as filtering, assembly, quality assessment, and graph visualization

1. Open SPAdes Assembly.
2. Select the required input files.
3. Choose the output directory.
4. Select the workflow steps and configure parameters.
5. Select Run Pipeline.
6. Record the Run ID shown by GenXFlow for later reference.

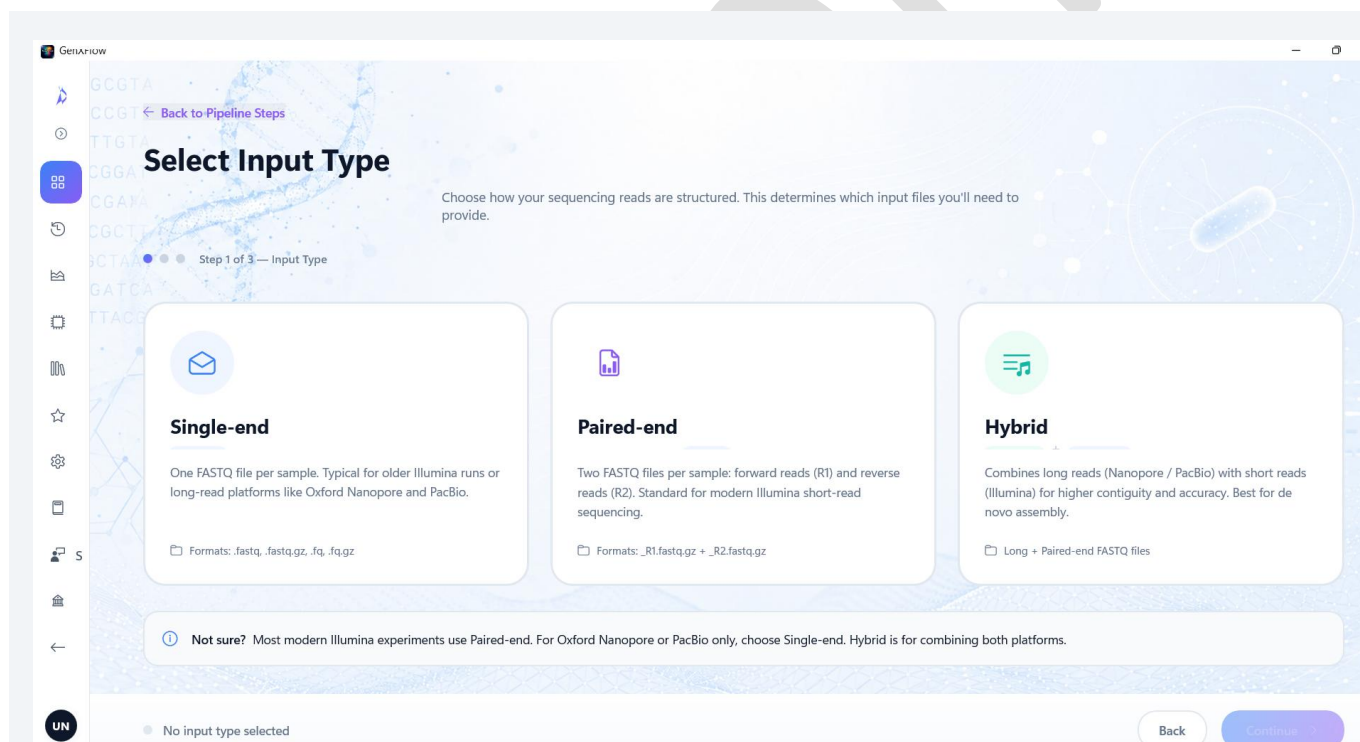


SCREENSHOT 9
SPAdes pipeline configuration

6.2 Flye Long-Read Assembly

Use Flye for long-read assembly workflows using Nanopore or PacBio data.

1. Open Flye Assembly.

2. Add one or more long-read FASTQ files.
3. Choose the output directory.
4. Select the pipeline steps.
5. Set genome size, minimum read length, retention percentage, and threads as required.
6. Select Run Pipeline.

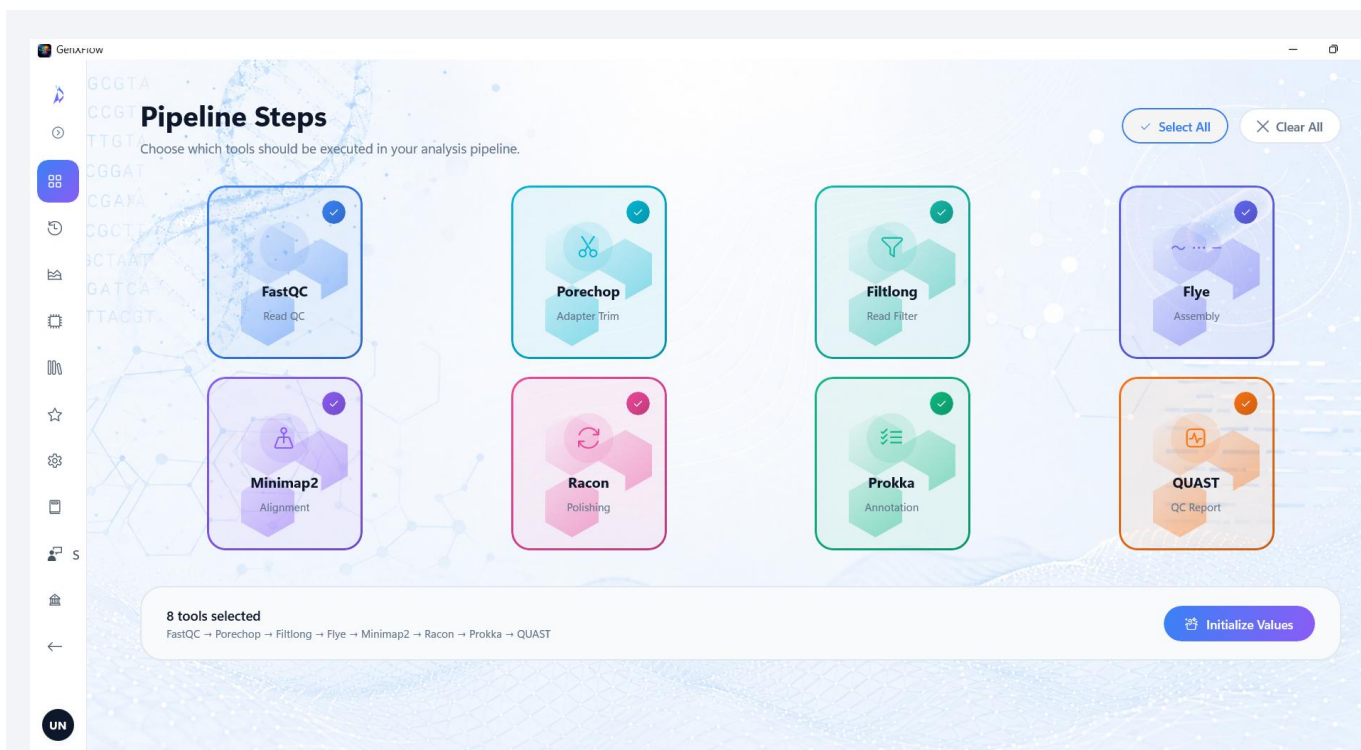

SCREENSHOT 10
Flye pipeline configuration


6.3 QIIME2 Amplicon / Metagenomics

Use QIIME2 for supported 16S, ITS, and related microbial community workflows.

Input / Parameter	Description
Input files	FASTQ files or QIIME2 artifacts (.qza), as supported
Metadata file	TSV/TXT sample metadata
Metadata column	Grouping column used for analysis
Target gene	16S, ITS, or 18S
Primers	Forward and reverse primer sequences where applicable
Truncation lengths	DADA2 read truncation settings
Sampling depth	Rarefaction depth for diversity analysis
Output directory	Location for results

1. Open QIIME2 Amplicon.

2. Add input data and the metadata file.
3. Select the metadata grouping column.
4. Configure primers and denoising/diversity parameters as required.
5. Select the required reference classifier/database.
6. Choose the output directory.
7. Select Run Pipeline.

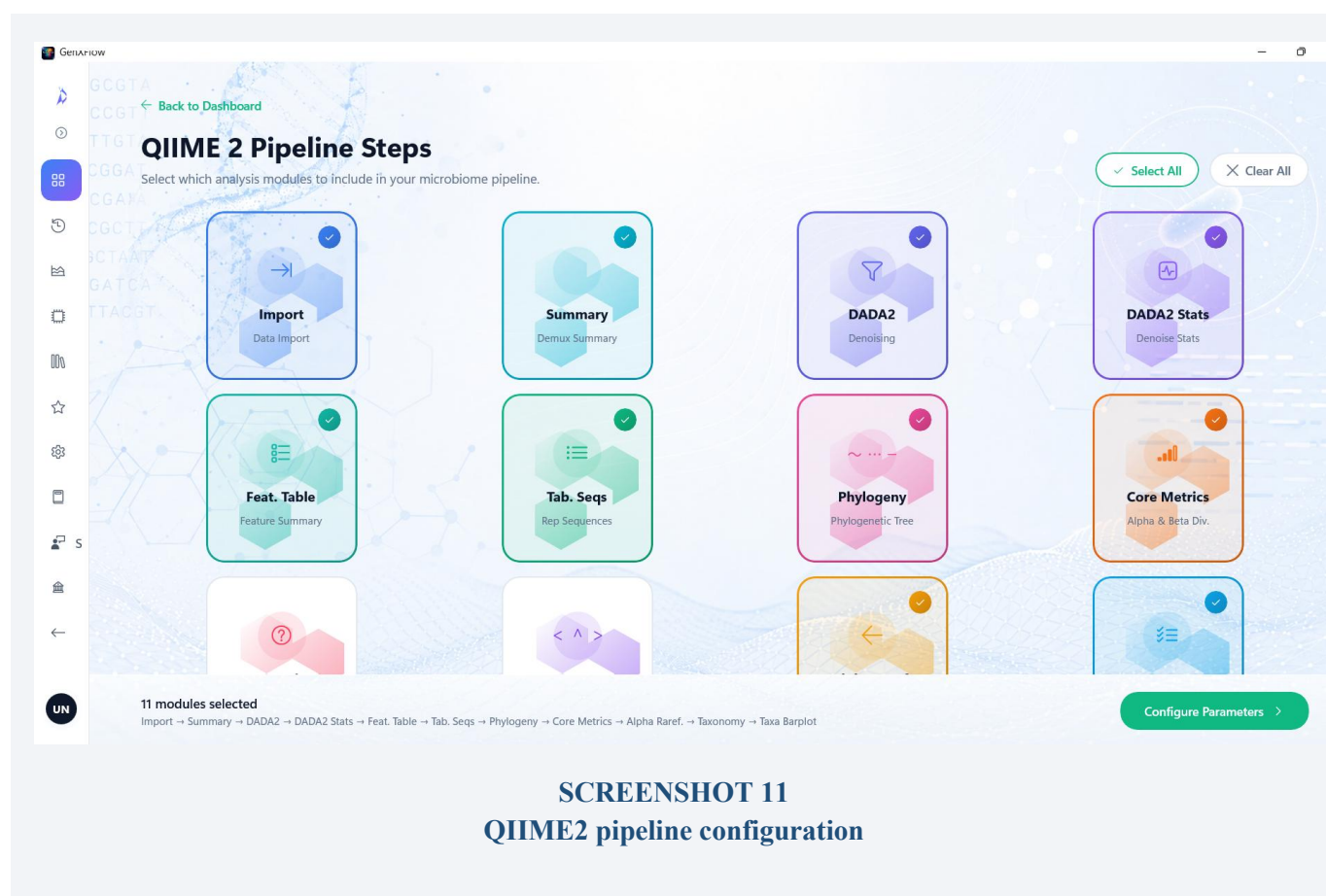


6.4 Sequence Similarity Analysis

Compares two FASTA files using exact matching, global alignment, or local alignment.

Setting	Options
FASTA File 1 / 2	Two FASTA inputs
Method	Exact, global, or local
Matching Mode	By sequence ID or by order
Output	Output directory and results file name

1. Open Similarity Analysis.
2. Select the two FASTA files.
3. Choose the comparison method and matching mode.
4. Set the output directory and result filename.

5. Select Run Analysis.

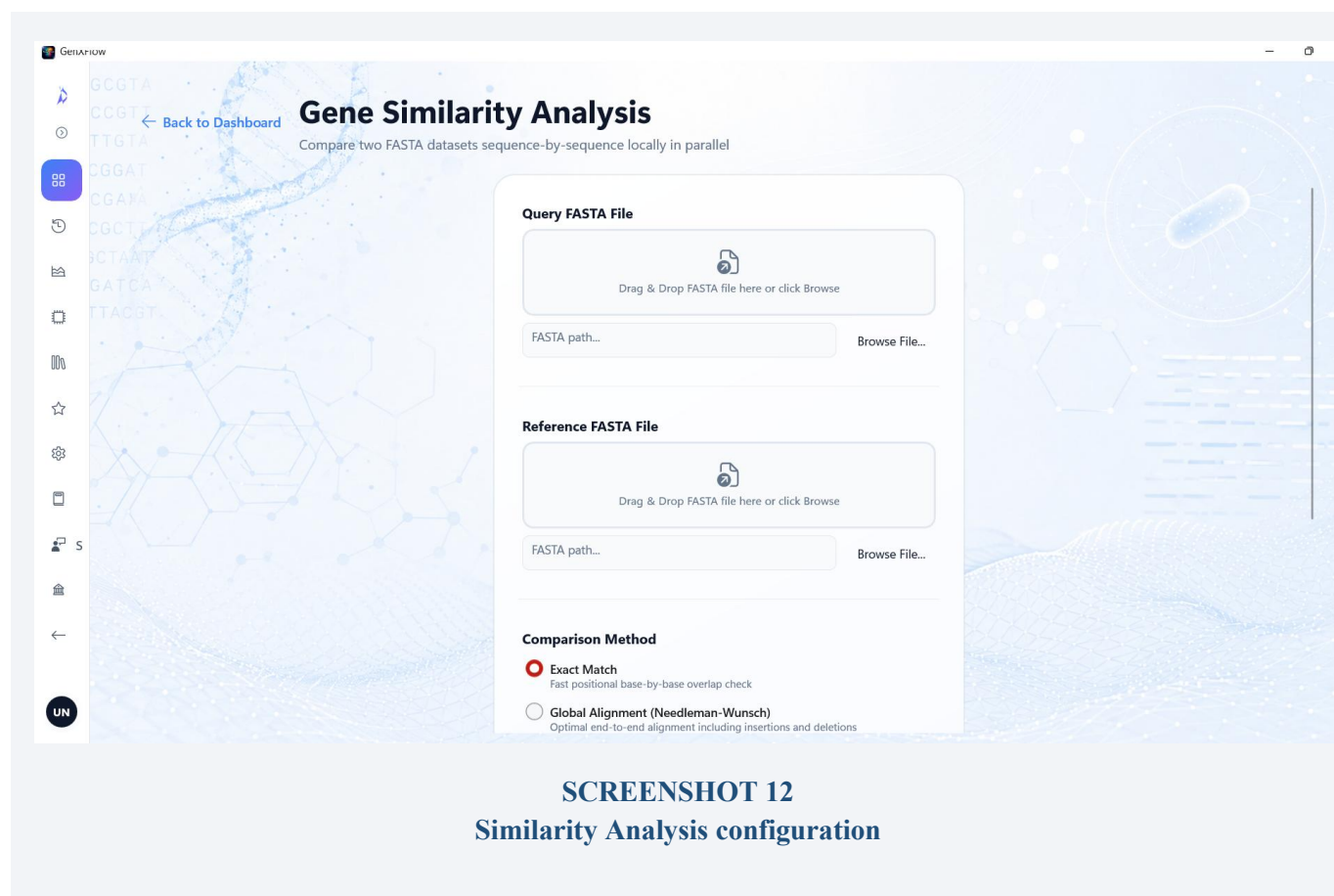

6.5 BLAST Search

Use BLAST Search to submit a nucleotide or protein sequence for remote homology searching through the supported NCBI BLAST service.

Setting	Typical value / description
Sequence	Raw sequence or FASTA-formatted sequence
Database	For example nt, nr, or another supported database
Program	blastn, blastp, blastx, tblastn, or tblastx as appropriate
E-value	Statistical significance threshold
Max Hits	Maximum number of results to retrieve
Email	Optional contact information for supported queries

1. Open BLAST Search.
2. Paste or enter the sequence.
3. Select the program and database.
4. Set the E-value and maximum hits.
5. Select Run BLAST.

GenXFlow

Dashboard

Runs

Analytics

Resource Monitor

Database Manager

License & Updates

Settings

User Guide

Support & Feedback

Activity Center

Notifications

universitypad.2023@g
Licensed Plan

Comparison Results

Run Folder: C:\ProgramData\GenXFlow\Runs\5b028d9-34ce-4fd4-ac94-353ce5f4bbfc

New Comparison

Matched Pairs: **102862**

Missing (Query/Ref): **0 / 740**

Average Similarity: **56.16%**

Highest Similarity: **100.00%**

Lowest Similarity: **19.11%**

Filters & Search

Search Seq ID

Query or Reference ID...

Min Similarity %: 0%

Min Identity %: 0%

☐ Has Mismatches

☐ Has Gaps

Export Options

Export Summary TSV

Export Summary CSV

#	Query ID	Reference ID	Length	Similarity	Identity	Gaps	Actions
1	SRR40321876.44312	SRR40321877.44312	301	55.48%	55.48%	0	View BLAST
2	SRR40321876.44312	SRR40321877.44312	293	43.69%	43.69%	0	View BLAST
3	SRR40321876.44458	SRR40321877.44458	301	56.15%	56.15%	0	View BLAST
4	SRR40321876.44458	SRR40321877.44458	293	31.06%	31.06%	0	View BLAST
5	SRR40321876.45292	SRR40321877.45292	301	60.13%	60.13%	0	View BLAST
6	SRR40321876.45292	SRR40321877.45292	293	42.32%	42.32%	0	View BLAST
7	SRR40321876.12757	SRR40321877.12757	293	45.73%	45.73%	0	View BLAST
8	SRR40321876.12757	SRR40321877.12757	301	64.45%	64.45%	0	View BLAST
9	SRR40321876.17800	SRR40321877.17800	293	43.69%	43.69%	0	View BLAST
10	SRR40321876.17800	SRR40321877.17800	301	92.36%	92.36%	0	View BLAST

Showing 102862 of 102862 matched pairs

SCREENSHOT 13
BLAST Search configuration

7. Monitoring & Run Management

7.1 Run progress

Pipeline Status
ID: cb321ebc-5efb-43bc-96f5-a0181e6cef95 Live

Analysis Running
Pipeline is actively executing...

Elapsed Time: **00:36** | Pipeline: **FLYE** | Log Lines: **10** | Status: **RUNNING**

Execution Progress
Step-by-step tool execution status

Tool	Status
porechop	Running
fllong	Pending
flye	Pending
minimap2	Pending
racon	Pending
fastqc	Pending
prokka	Pending
quast	Pending

Execution Log (Auto-scroll ON, Copy)


```

290 / 1,956 (14.8%)
300 / 1,956 (15.3%)
310 / 1,956 (15.8%)
320 / 1,956 (16.4%)
330 / 1,956 (16.9%)
340 / 1,956 (17.4%)
350 / 1,956 (17.9%)
360 / 1,956 (18.4%)
370 / 1,956 (18.9%)
380 / 1,956 (19.4%)
390 / 1,956 (19.9%)
400 / 1,956 (20.4%)
410 / 1,956 (21.0%)
  
```

SCREENSHOT 14
Run progress / monitor page

The monitor shows the status of the run, the current workflow step, completed steps, progress percentage, elapsed time, and available execution information.

7.2 Run History

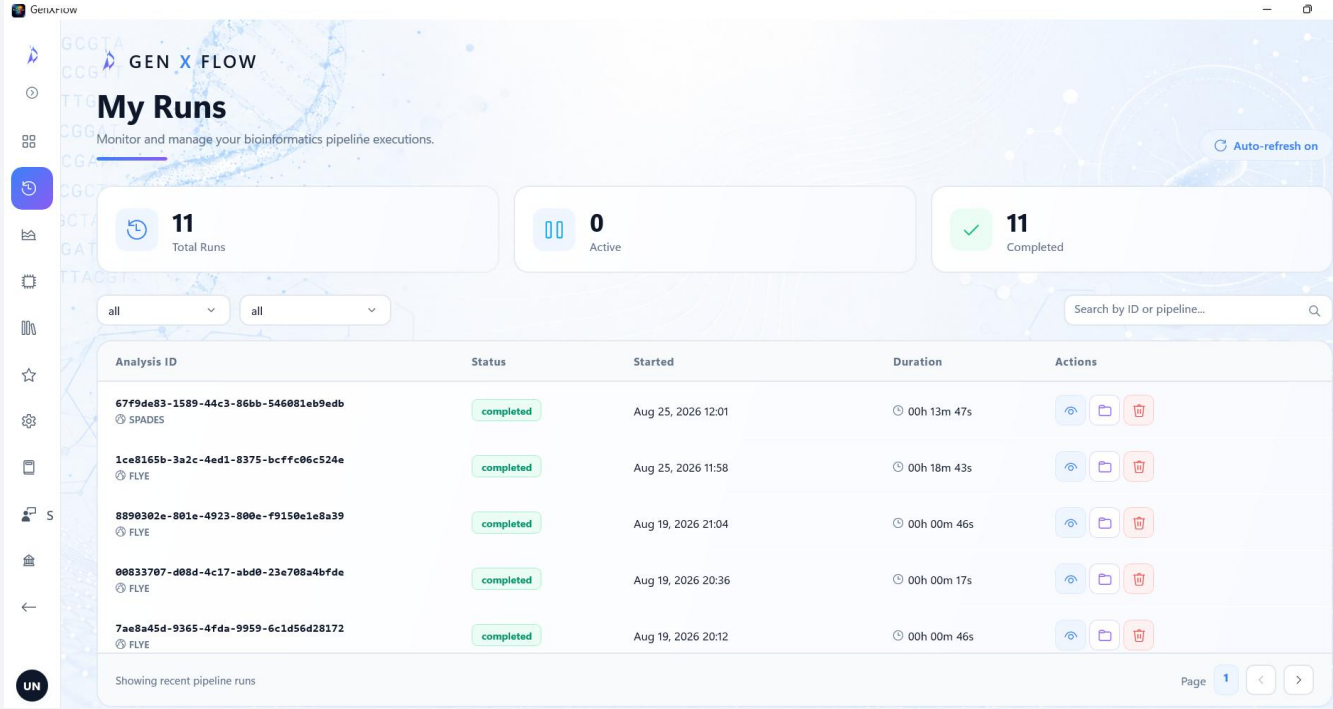

The screenshot displays the 'My Runs' interface in GenXFlow. At the top, it says 'GEN X FLOW' and 'My Runs' with a subtitle 'Monitor and manage your bioinformatics pipeline executions.' There are three summary cards: '11 Total Runs', '0 Active', and '11 Completed'. Below these are two dropdown menus set to 'all' and a search bar labeled 'Search by ID or pipeline...'. The main part of the page is a table with the following data:


Analysis ID	Status	Started	Duration	Actions
67f9de83-1589-44c3-86bb-546081eb9edb SPADES	completed	Aug 25, 2026 12:01	00h 13m 47s	[View] [Download] [Delete]
1ce8165b-3a2c-4ed1-8375-bcffc06c524e FLYE	completed	Aug 25, 2026 11:58	00h 18m 43s	[View] [Download] [Delete]
8890302e-801e-4923-800e-f9150e1e8a39 FLYE	completed	Aug 19, 2026 21:04	00h 00m 46s	[View] [Download] [Delete]
00833707-d08d-4c17-abd0-23e708a4bfde FLYE	completed	Aug 19, 2026 20:36	00h 00m 17s	[View] [Download] [Delete]
7ae8a45d-9365-4fda-9959-6c1d56d28172 FLYE	completed	Aug 19, 2026 20:12	00h 00m 46s	[View] [Download] [Delete]

At the bottom of the table, it says 'Showing recent pipeline runs' and 'Page 1' with navigation arrows.

SCREENSHOT 15
Run History page

Use Run History to find current and previous executions. Search by Run ID and filter by status or pipeline type.

7.3 Run details

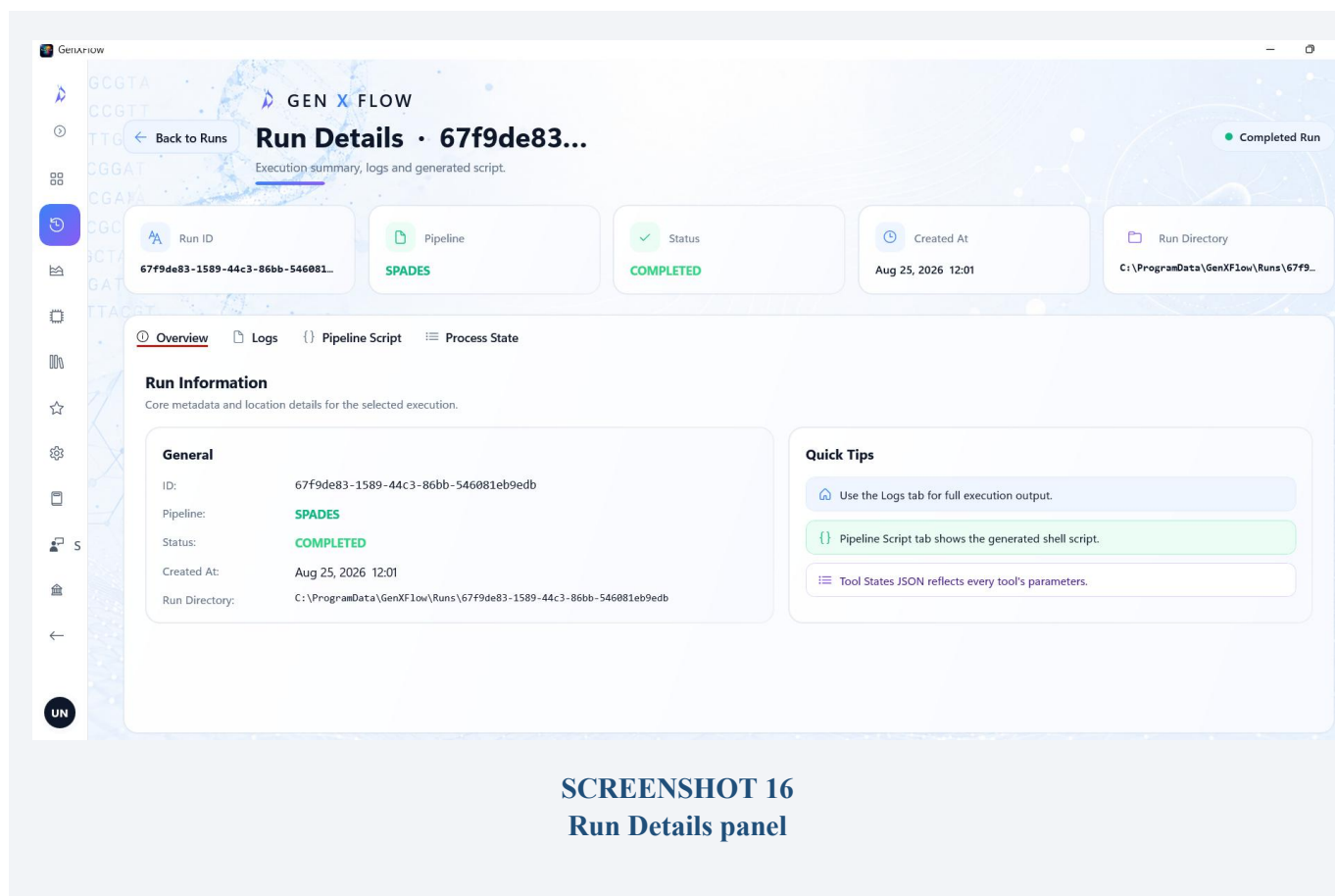

SCREENSHOT 16
Run Details panel

7.4 Cancelling a run

1. Open Run History or the active run monitor.
2. Locate the running pipeline.
3. Select Cancel.
4. Confirm the cancellation if prompted.


Cancellation may take some time to take effect because the active workflow step may need to finish before the run stops.

7.5 Deleting a run

1. Open Run History.
2. Select the run to remove.
3. Select Delete and confirm the action.

Deleting a run removes the run record from GenXFlow. Output files may remain on disk and should be removed separately when no longer required.

8. Reference Databases

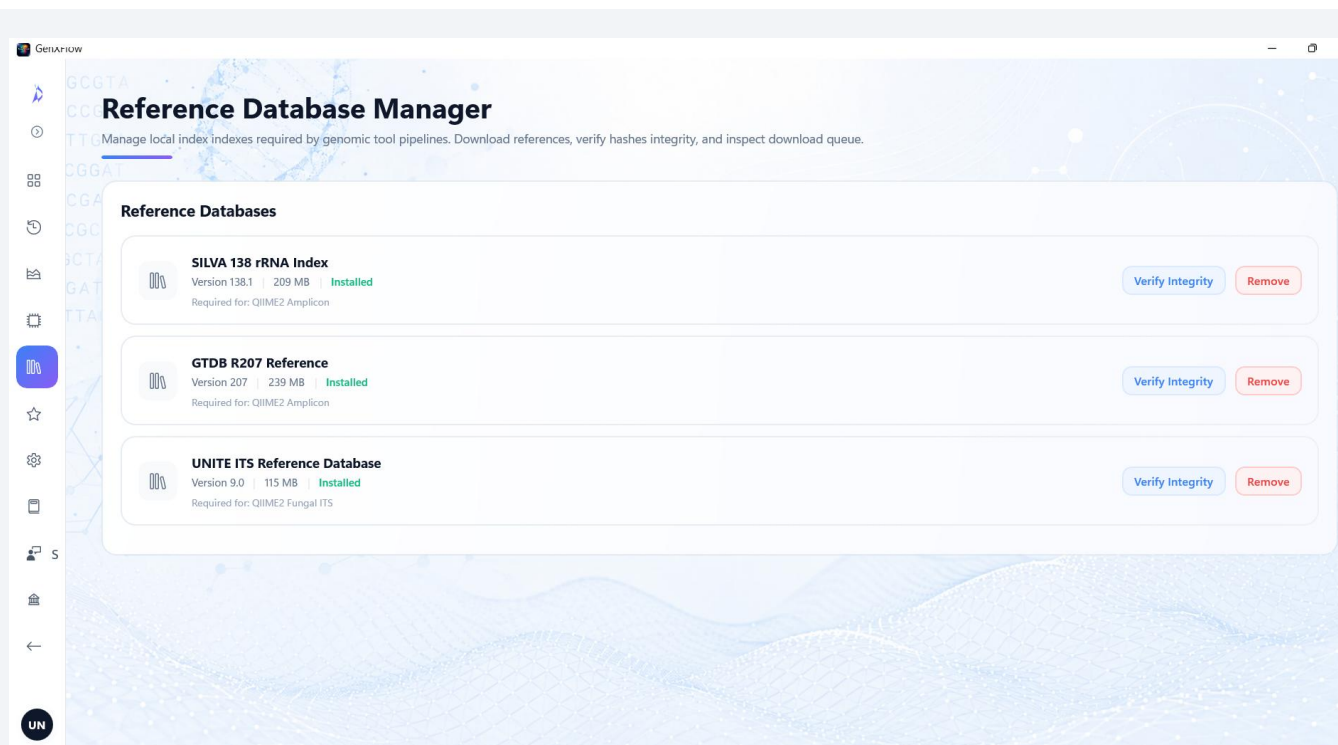


SCREENSHOT 17
Database Manager page

The Database Manager is used to download and verify reference databases needed by supported workflows, particularly QIIME2 taxonomy classification.

Database	Use
SILVA 138 rRNA Index	QIIME2 16S classification
GTDB R207 Reference	Supported taxonomy / profiling workflows
UNITE ITS Reference	QIIME2 fungal ITS classification

8.1 Download a database

1. Open Databases.
2. Locate the required database.
3. Select Download.
4. Wait for the download and integrity verification to complete.
5. Confirm the status changes to Installed.

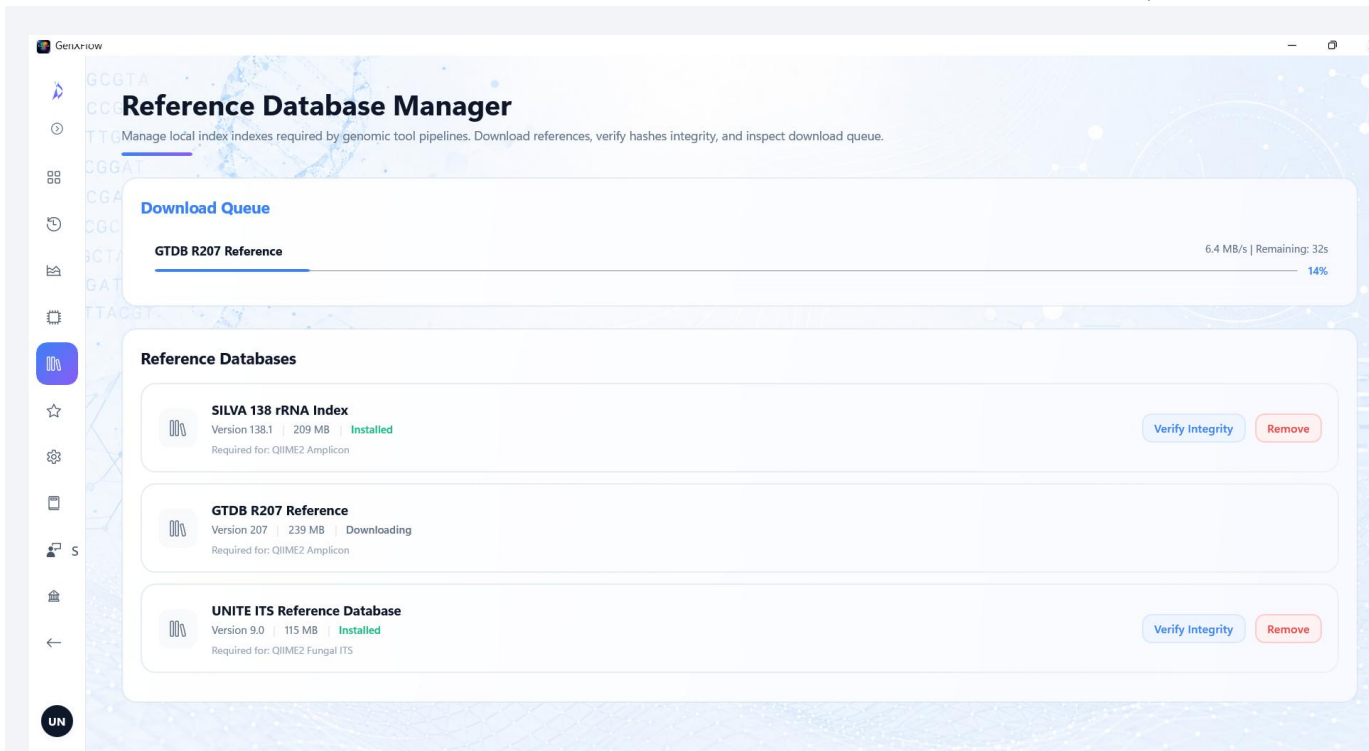


SCREENSHOT 18
Database download progress

8.2 Verify or remove a database

- Use Verify to check the integrity of an installed database.
- Use Remove to delete a database that is no longer needed.
- If verification fails, re-download the database.

9. Notifications & System Health


9.1 Notifications

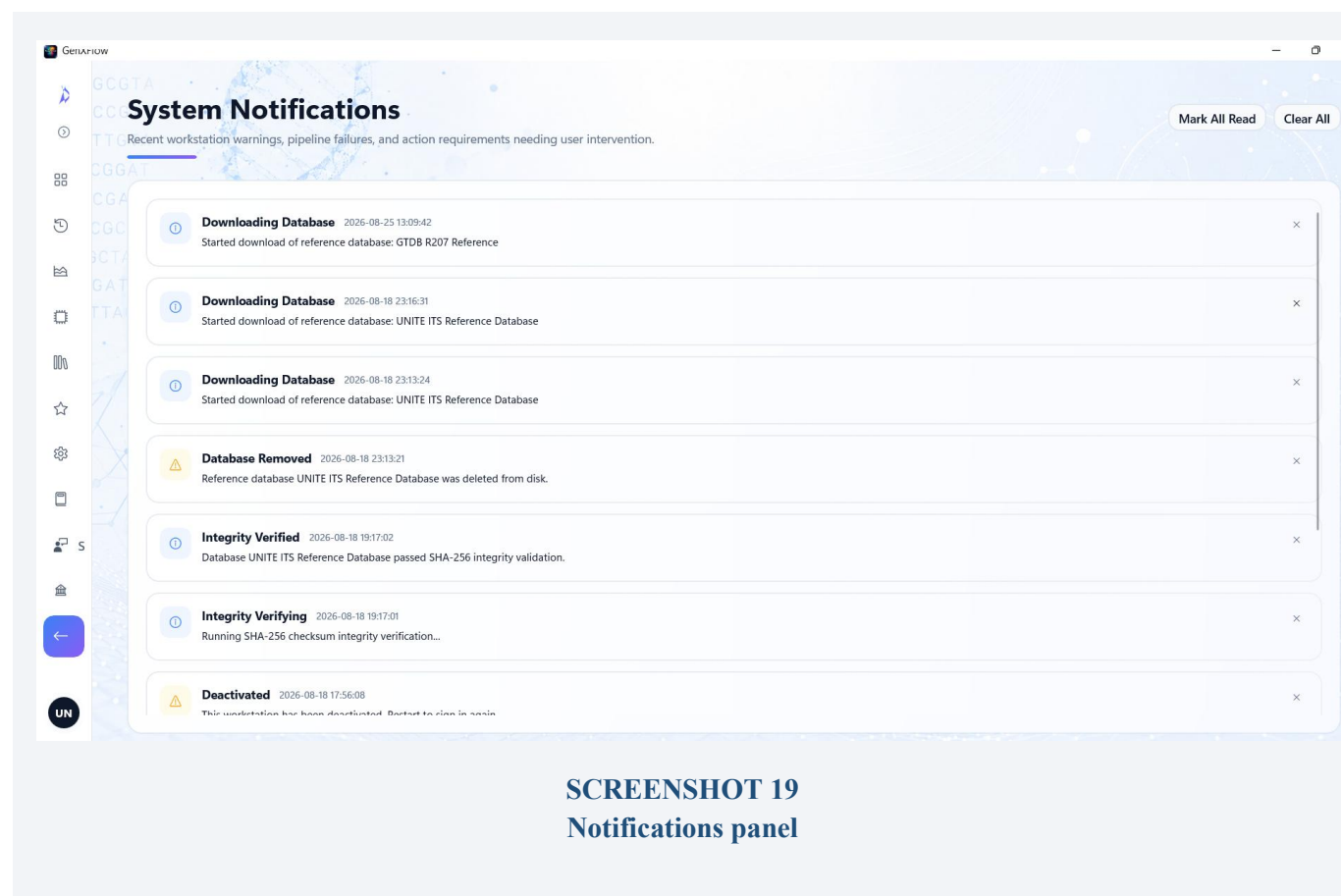

SCREENSHOT 19
Notifications panel


The notification center reports important application, pipeline, licensing, and update events. Use Mark as Read, Mark All as Read, Dismiss, or Clear All as appropriate.

9.2 System health

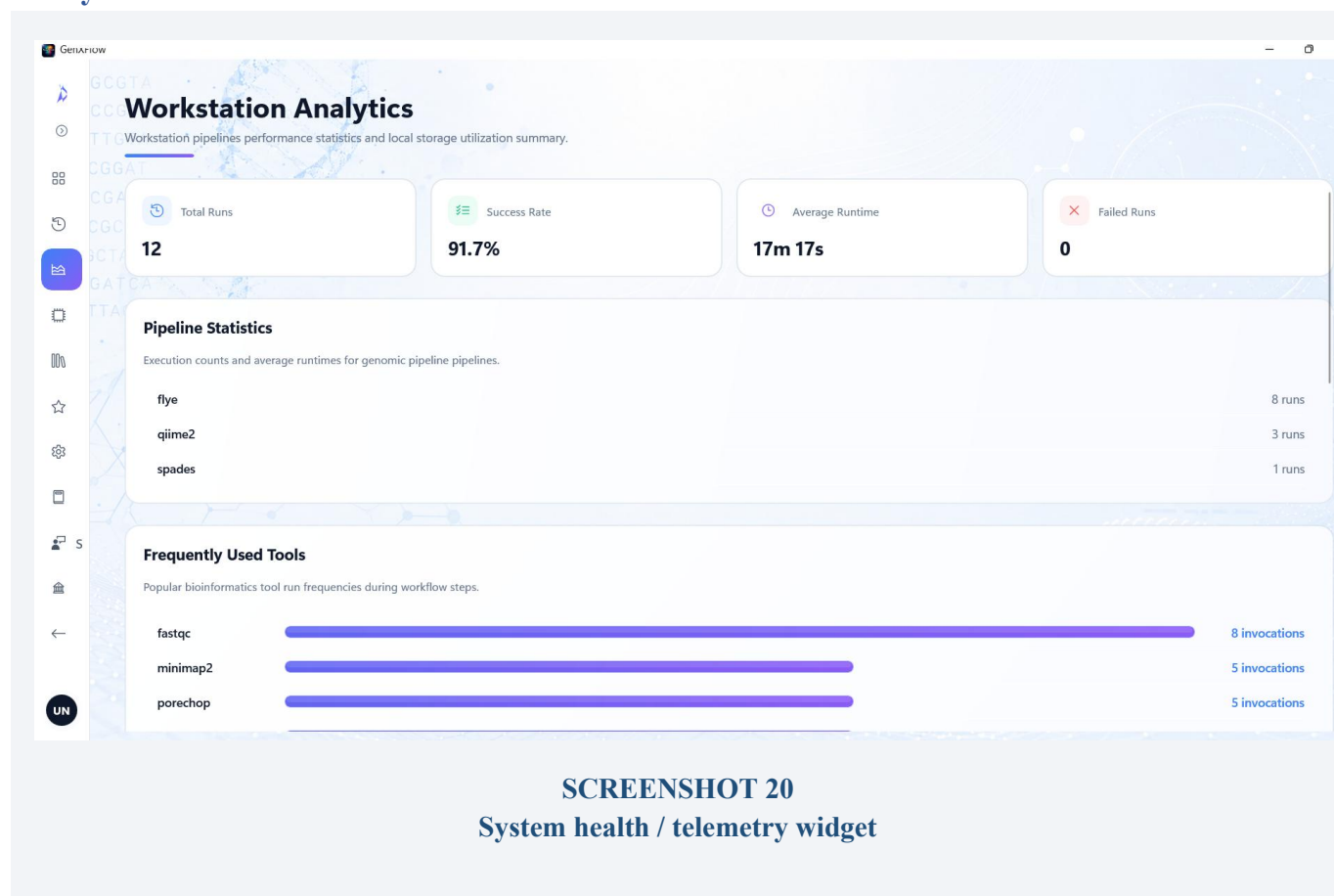

Indicator	Customer action
RAM available is low	Close unnecessary applications before large analyses
Free disk space is low	Free or allocate additional storage before running large workflows
Docker is not running	Start or restart Docker Desktop
CPU usage is very high	Allow current workloads to finish before starting another large run

10. Application Updates

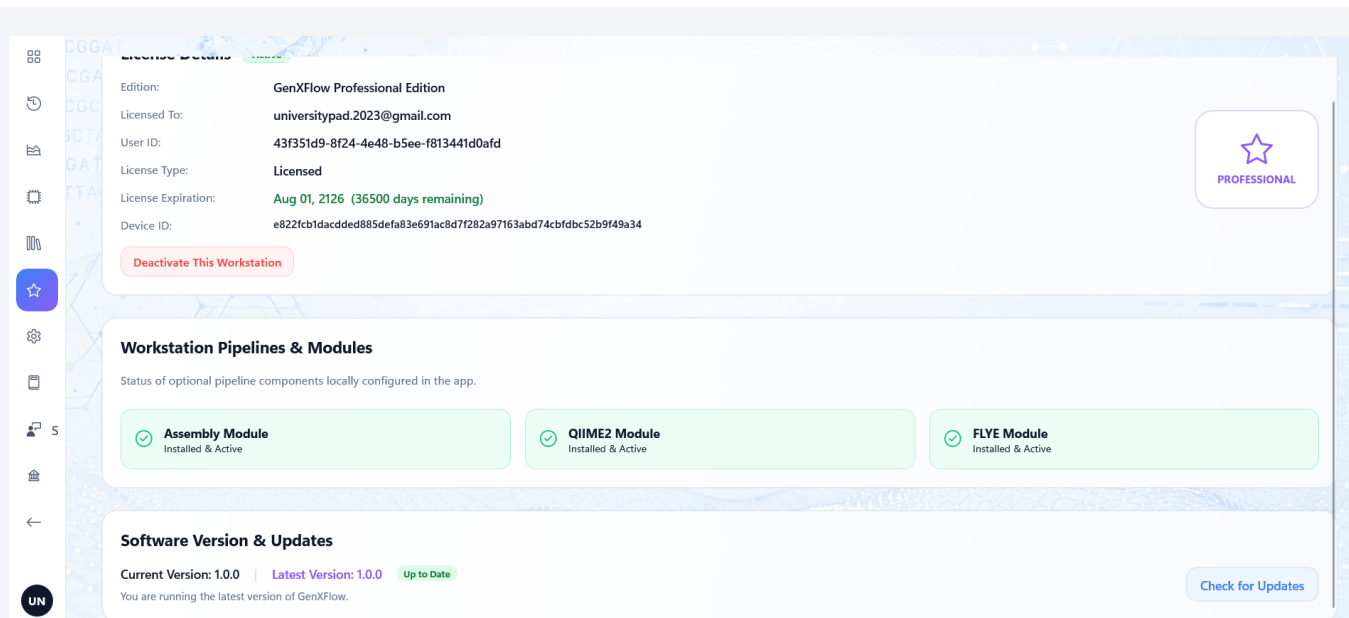

SCREENSHOT 21
Update notification / update banner

GenXFlow checks for available application updates. When an update is available, follow the prompt to download and install the new version.


1. Review the update notification.
2. Select Download Update.
3. Wait for the download and verification to complete.
4. Select Install Now when prompted.
5. Allow GenXFlow to close and relaunch as part of the update process.

11. Feedback & Troubleshooting

11.1 Submit feedback

Use Settings → Feedback, or the Feedback control provided in your GenXFlow build, to report a bug, request a feature, or send general feedback.

1. Select the feedback type.
2. Describe the issue or request clearly.
3. Add a contact email if follow-up is needed.
4. Select Submit.

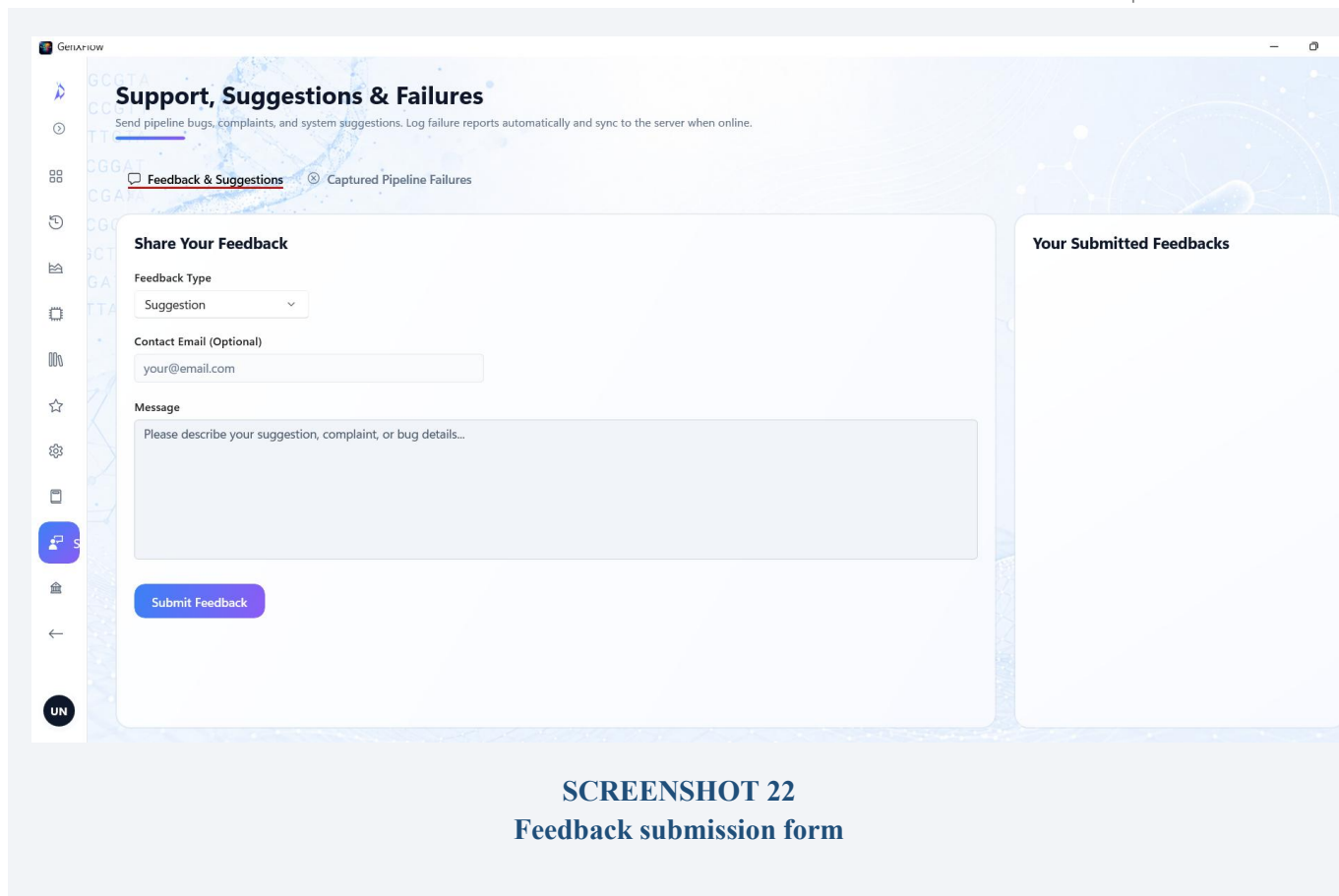

SCREENSHOT 22
Feedback submission form

11.2 Common issues

Issue	What to check
Pipeline will not start	Sign in and verify your license; make sure required inputs are present; confirm Docker is available for Docker-based workflows.
Docker-related error	Start Docker Desktop and retry. If the issue continues, restart Docker Desktop.
Trial limit reached	A paid license is required for additional pipeline runs.
Run remains at 0%	Check Docker availability, internet access, and whether required images/resources are available.
Database download fails	Check network connectivity and retry the download.
Database verification fails	Remove or cancel the incomplete database and download it again.
Cannot sign in	Check internet connectivity and verify your account credentials.
QIIME2 taxonomy fails	Confirm that the required SILVA or UNITE database has been downloaded and installed.
Insufficient disk space	Free disk space or use an appropriate storage location for outputs and databases.

12. Support Notes

When contacting ArqGene support, provide the following information where possible:

- GenXFlow version.
- Operating system version.
- Pipeline name.
- Run ID.
- A short description of the issue and the step at which it occurred.
- Relevant error message or screenshot.

Use the in-app Feedback function or the support channel provided with your customer license for assistance.

Appendix A. Screenshot Checklist

Paste actual product screenshots into the numbered locations above. Capture the full application area where practical, and make sure customer-sensitive information such as personal email addresses, passwords, API keys, or private sample identifiers is removed or blurred.

#	Screenshot	What to capture
1	Product cover / welcome	Optional branded product image.
2	Installation wizard	Installer welcome/first step.
3	First launch	Loading or initialization screen.
4	Login	Email/password sign-in page.
5	Registration	Create Account form.
6	Forgot Password	Password recovery form.
7	License information	Status, edition, expiry and account details.
8	Main dashboard	Pipeline cards, run stats, recent runs, health, notifications.
9	SPAdes configuration	Inputs, tools, parameters, output directory, Run button.
10	Flye configuration	Long-read inputs, parameters, output directory, Run button.
11	QIIME2 configuration	Inputs, metadata, parameters, database and Run button.
12	Similarity Analysis	FASTA inputs, method, matching mode, output.
13	BLAST Search	Sequence, program, database, thresholds, Run button.
14	Run monitor	Status, current step, progress and log.
15	Run History	Run list, filters and statuses.
16	Run Details	Selected run details and accessible logs/output information.
17	Database Manager	Database cards and statuses.
18	Database download	Active download and verification progress.
19	Notifications	Notification center and examples.
20	System health	CPU, RAM, disk and Docker status.
21	Update notification	New version notification and download action.
22	Feedback form	Feedback category, message and submit action.

End of Customer User Manual